High Power Fractional CO₂ Laser plus Fractional RF Micro-needling

The Advancement of **NEW** Laser Technology

RF Micro-needling

FRA EIS

FRAXIS DUO is a combination CO₂ laser and microneedle fractional RF technology for the treatment of scars and photo-aging. CO₂ fractionals enables intense superficial stimulation and micro-needle fractional RF provides higher volumetric heating and deeper heat diffusion.

Face

- Skin Imperfections
- Pore Reduction
- Wrinkle Reduction
- Skin Restoration
- Acne Scar
- Resurfacing

Body

- Scars
- Stretch Marks

Women's

Health

CO₂ Laser

Keratosis Pilaris

Women's Health

- Vaginal Health
- Collagen Growth







FRA:S

High Power Fractional CO2 Laser Plus RF Micro-needling

> Multiple applications, combination therapy

The benefits of combining a fractional CO₂ laser RF micro-needling provides the patient with both superficial and deeper treatments. This multi treatment approach is beneficial for many different skin abnormalities and treats a larger range of applications.

Maximum efficacy and safety

With dual energies and various modes, it enables medical professionals to easily select the most effective mode and energy. This allows the practitioner the ability to treat a variety of conditions, safely targeting specific skin layers without damage to the surrounding tissue for patient comfort and safety.

Cost effective and versatile


The combination modality provides maximum treatment flexibility for photo-aging, scars to striae distensae. Such a dual multi-functional system makes it easy to expand patient services without the need to purchase additional equipment.

Smart and easy to use

User friendly designed 10.4" full color touchscreen helps physicians easily choose the optimal parameter and a weight type articulated arm provides accurate movement in any position to easily place the beam centering and minimize physician fatigue.

Combination therapy for optimal treatment result

Combining 2 techniques optimize the result for certain indications that require both superficial and deeper heating coagulation.

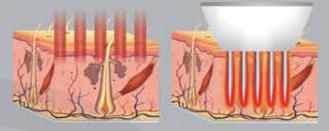

CO² laser mode

Micro-needle RF mode

CO₂ Laser Handpiece





Non-surgical laser treatment for vaginal health.

▲ Disposable external slide



SPECIFICATIONS

CO₂ Fractional Laser System

Wave length	10,600nm (TEM00)	
Laser power	Up to 30W (Up to 150mJ)	
Pulse duration	20 - 5,000us	
Overlap	1 - 10th	
Pixel size	≥100 micron	
Distance	0.1 - 2.0mm	
Treatment area	1cm x 1cm - 20 x 20mm	
Pixel quantity	Up to 40, 401	
Repetition	0.2 - 2 Sec/Single	
Surgical	0.2/0.3/0.5/0.7/	
handpiece (option)	0.9/1.1/1.3mm	

Micro-needle Fractional RF System

RF Frequency	2MHz± 10%	INVISIBLE AND/OR VISIBLE LASER RADAITION WORK DE YOR SKILL RESOLUCE TO RIBECT OF SCATTERED RADAITION MUNICIPALITY STATEMENT MUNICIPALITY STATEMENT MUNICIPALITY MUNICIPALITY STATEMENT MUNICIPALITY MUNICIPALITY MUNICIPALITY MUNICIPALITY
RE Power	Up to 25W±10%	
Treatment duration (RF emission)	50ms ~950ms	CLASS IV LASER PRODUCT
Intensity	0-100 Level (2 / 5 / 10 Step)	
Needle thickness	0.25mm	
Needle figure	25 electrode needle (5x5)	
	10 electrode needle (1x10) (Option)	
Treatment area	25Pin (7.8 x 7.8mm) / 10Pin (0.25 x 16mm)	
Adjustable depth	0.5 ~ 3.5mm (0.1 Step)	
Repetition	0.2 / 0.5 / 0.8 / 1 / 2 Sec, Sir	ngle
Electrical power	100 - 240VAC, 50 / 60Hz	



LIVE DEMOS ON DEMAND! What are you waiting for? (888) 527-7368 | office@laserent.net

DANGER